

Q-NADMED BLOOD

Quantitative assay kit for NAD⁺ and NADH in Blood

FOR SINGLE USE ONLY

These instructions must be read in its entirety before using this product.

CE  FOR *IN VITRO* DIAGNOSTIC USE

GENERAL INFORMATION

- A. Proprietary name: Q-NADMED Blood: quantitative assay kit for NAD⁺ and NADH in Blood
- B. Catalog number: IVD_001, 40 samples (96-well format)
- C. IFU issued: May 2022

Manufacturer:

NADMED Ltd / Oy

Haartmaninkatu 4

00290 Helsinki, FINLAND

www.nadmed.fi, info@nadmed.fi

TABLE OF CONTENTS

GENERAL INFORMATION	2
TABLE OF CONTENTS	3
INTENDED USE	4
CLINICAL BACKGROUND	4
PRINCIPLES OF THE ASSAY	4
REAGENTS PROVIDED	5
PRECAUTIONS AND WARNINGS	5
STORAGE AND EXPIRATION DATING OF REAGENTS	6
OTHER MATERIAL REQUIRED	6
REAGENT PREPARATION	7
SAMPLE COLLECTION AND STORAGE	7
PRACTICAL CONSIDERATIONS	7
EXTRACTION OF NAD⁺ AND NADH FROM BLOOD	9
PREPARATION OF STANDARDS	10
ASSAY PROCEDURE	11
CALCULATION OF RESULTS	12
INTERNAL QUALITY CONTROL	12
PERFORMANCE AND LIMITATIONS	13
SYMBOLS	17
SCHEMATIC PICTURES	19
NOTES	20
PLATE LAYOUT	21

INTENDED USE

Q-NADMED Blood, an *in-vitro* diagnostic medical device, is an analytical assay kit for measurement of concentrations of NAD⁺ and NADH metabolites in human whole blood. The assay is quantitative. The intended users of the Q-NADMED assay kit are trained laboratory personnel. The first intended purpose is to detect systemic changes of NAD⁺ and NADH. The primary intended users of the assay results are healthcare professionals who interpret the obtained results in the context of a disease/health status. Results of the Q-NADMED assay kit can be used for decision making on treatment such as supplementation with NAD precursors. The second intended purpose of Q-NADMED assay kit is to monitor NAD⁺ and NADH levels in patients receiving treatment such as supplementation of NAD precursors and to adjust the dose.

CLINICAL BACKGROUND

NAD⁺ and NADH metabolites are involved in adaptation of human body metabolism and energy homeostasis to changing endogenous and exogenous conditions. Accumulated research data show that systemic levels of NAD⁺ decrease in response to manifested disease, signaling of an imbalance of body energy homeostasis. Degree of NAD⁺ decrease varies in different patients and in different pathologies. Progressive decline of NAD⁺ levels makes it impossible for the body to maintain its basic metabolic functions to survive even in conditions of ongoing therapy. Q-NADMED allows screening patients for NAD⁺ and NADH deficiency with the aim to correct it and increase efficiency of treatment. Ongoing research on the contribution of NAD⁺ and NADH to mechanisms and progression of different diseases is very active. List of pathologies with suspected changes of NAD⁺ and NADH concentrations is constantly expanding with already published evidence for mitochondrial disease, aging, sepsis, viral infections, cardiovascular and kidney disease, Diabetes type I and II, neurological disorders and cancer.

PRINCIPLES OF THE ASSAY

The principle of the assay is a cyclic enzymatic reaction with a colorimetric end-point detection. First, NAD⁺ and NADH metabolites are extracted together from a blood sample in a single step. Then, the extract is divided into two parts. In the first part, NAD⁺ is stabilized while NADH is removed, whereas in the second part, NADH is stabilized while NAD⁺ is removed. Next, NAD⁺ and NADH metabolites are analyzed on two separate plates by an enzymatic reaction coupled to a color change. The intensity of the color change in the assay is linearly proportional to the concentration of NAD⁺ or NADH in the reaction mixture.

REAGENTS PROVIDED

REAGENTS	DESCRIPTION	BEFORE THE ASSAY
BUFFER A	28 mL of Extraction Buffer	Equilibrate to room temperature (15-25°C).
NAD⁺ stabilizing reagent	8 mL of buffered solution for NAD ⁺ measurement	Equilibrate to room temperature (15-25°C).
NADH stabilizing reagent	8 mL of buffered solution for NADH measurement	Equilibrate to room temperature (15-25°C).
NAD⁺ standard stock	1 vial with 40 µL of 1 mM NAD ⁺	See preparation guide.
NADH standard stock	1 vial with 40 µL of 1 mM NADH	See preparation guide.
BUFFER C	2 vials with 19 mL of Assay Buffer	Equilibrate to room temperature (15-25°C).
Assay color reagent	2 vials with 3 mL of reagent with the assay color.	Equilibrate to room temperature (15-25°C). Should be used within 1 h after equilibration.
Enzyme	2 vials with 40 µL of enzyme. One per plate.	Thaw only before adding into Master Mix.
Stop Solution	3 mL of solution to stop the assay reaction.	Equilibrate to room temperature (15-25°C).
Positive control	1 vial with 150 µL of NAD ⁺ /NADH solution	Thaw on ice-water bath for 10 min prior to extraction

PRECAUTIONS AND WARNINGS

Safety

For *in vitro* diagnostic use only by trained personnel.

The Stop Solution may cause skin, eye and respiratory irritation. Avoid breathing in fumes.

Assay color reagent may cause skin irritation. Handle with care, use gloves.

BUFFER A can cause eye irritation. Handle with care, use goggles.

Do not smoke, drink, eat or apply cosmetics in the working area. Wear protective gloves, clothing, and eye protection. Wash hands thoroughly after handling.

Identified hazards of the chemicals presented in this kit and the appropriate warning information associated with those hazards are presented in Q-NADMED Safety Data Sheet (SDS).

Disposal of used kit components is described in in Q-NADMED Safety Data Sheet (SDS).

STORAGE AND EXPIRATION DATING OF REAGENTS

- Before opening, all kit components are stable until the expiry date, indicated on the label, if kept at -20°C. If stored at -80°C, the components are stable for 6 months. Avoid temperature fluctuations in the freezer.
- After thawing, BUFFER A, BUFFER C, NAD⁺ and NADH Stabilizing Reagents and Stop Solution are stable for 2 weeks at room temperature
- Assay Color reagent is stable for 1h at room temperature after thawing
- The Enzyme should be used directly after thawing
- Standards and the Positive Control should be prepared and used on the day of the assay
- Standards should be protected from light

OTHER MATERIAL REQUIRED

The following materials are required but not provided in the kit:

- MilliQ water
- Test tubes (1.5 mL, e.g., microcentrifuge tubes) for dilution of standards and sample preparation
- Two Multichannel Pipette reservoirs, one for pipetting Master mix, the other for pipetting Stop Solution
- Two 15 mL tubes (e.g., Falcons) to prepare the Buffer NAD⁺ and Buffer NADH for dilution of Standard stocks
- Two 96-well transparent polystyrene microplates with medium protein binding affinity suitable for colorimetric and absorbance assays.
- Heat block with adjustable temperature (up to 80°C)
- Table-top cooling microcentrifuge (max speed 20 000 x g)
- Ice
- Single channel pipettes (0.5-10 µL, 5-50 µL, 20-200 µL, 100-1000 µL) and multichannel pipettes (5-50 µL, 30-300 µL) and pipette tips, low retention
- Microplate reader capable of measuring absorbance at 570-573 nm
- Aluminum foil to protect tubes and plates from light

REAGENT PREPARATION

Bring all components for NAD⁺ and/or NADH to room temperature before use and keep at the temperature indicated in the table above before the assay.

If a precipitate forms in the buffers during storage, it should be redissolved by incubating the buffers at 37°C, and then cooled back down to 25°C before use. Do not shake Stop Solution vigorously.

- A. **Buffer A** - Ready to use
- B. **NAD⁺ stabilizing reagent** - Ready to use
- C. **NADH stabilizing reagent** - Ready to use
- D. **Enzyme** - Ready to use
- E. **Stop Solution** - Ready to use
- F. **50 µM NAD⁺ standard working stock**- adding 25 µL of 1 mM NAD⁺ standard stock (provided) into 475 µL of MilliQ water, vortex. Protect from light. (For further instructions see preparation of standards)
- G. **10 µM NADH standard working stock**- adding 10 µL of 1 mM NADH standard stock (provided) into 990 µL of MilliQ water, vortex. Protect from light. (For further instructions see preparation of standards)
- H. **Master Mix** – 1 vial of **Assay Color Reagent** (3 mL) should be mixed into 1 vial of **BUFFER C** to create 1 bottle of Master Mix required for 1 plate of either NAD⁺ or NADH assay. Protect from light by keeping it in the original amber bottle. Do not shake vigorously.

SAMPLE COLLECTION AND STORAGE

The kit is only suitable for whole blood samples. Whole blood samples should be collected into either heparin or EDTA collection tubes (e.g. BD K2E Vacutainer) and properly mixed. Samples can be analyzed fresh or frozen. Fresh blood can be stored for 3 hours at room temperature or at 4°C for up to 24 hours before analysis. Fresh samples can be aliquoted (e.g. 120 µL) and stored frozen at -20°C for two weeks or at -80°C for approximately a year. Frozen samples must be kept frozen at all times before the assay. Time intervals during sample handling should be consistent for all samples. For measurement of NAD⁺ and NADH, 100 µL of the whole blood is needed. Subsequent freeze-thaw cycles are not allowed.

PRACTICAL CONSIDERATIONS

- Do not use kit components beyond the expiry date.
- Do not mix materials from different kit lots. Subsequent freeze-thaw cycles are not allowed.
- This assay is NOT suitable for measurement of NAD⁺ and NADH in plasma or serum, cultured cells or tissues (human or animal).
- Thoroughly mix all reagents by gentle swirling. A quick centrifugation of smaller vials prior to opening is recommended.

- Both fresh and frozen blood samples can be used for the assay.
- We recommend extracting a maximum of eight samples at a time to minimize handling time.
- The analyses of NAD⁺ and NADH are done on two separate plates. We recommend performing the NAD⁺ and NADH assays on the day of extraction.
- To avoid cross-contamination, change to new pipette tips between the additions of each standard, sample, and reagent. Also, use separate reservoirs for Master Mix and Stop Solution.
- High precision pipettes and beveled tips with less retention will improve the precision.
- Assay Color Reagent is a yellow, light sensitive compound which turns brown upon enzymatic reaction. The exposure to direct sunlight or artificial light causes its unspecific color change to green. To minimize the light interference with the assay, we recommend
 - pipetting of Master Mix and then Stop Solution into the plate in the conditions of natural indirect light. To achieve these conditions, switch off artificial indoor lights and avoid pipetting close to the window.
 - covering the 96-well plate with a lid from an aluminum foil during the reaction time after addition of Master mix, as well as during the transfer to the plate reader after addition of Stop Solution. Do not wrap the plate with the aluminum foil due to risk of contamination and sample mixing during unwrapping.
- Positive control is an artificial sample containing known amounts of pure NAD⁺ and NADH.
- Technical requirements for spectrophotometric plate reader:
 1. measurement of light absorption at 570-573 nm,
 2. Option to adjust the scanning light brightness/intensity to **low**. In some plate readers, the brightness can be adjusted as the number of flashes per one measurement. In the latter case, set the number of flashes in the range from 5 to 10.
- Master Mix and Stop Solution contain detergents. To avoid bubbles, pipette Master Mix and Stop Solution by pressing the pipette to the first stop position, remove any bubbles in the wells with a small needle. Avoid touching the content of the wells by pipette tips.
- We recommend the following order of steps:
 1. Bring BUFFER A, NAD⁺ Stabilizing Reagent, NADH Stabilizing Reagent and Stop Solution to room temperature the day before the assay. These Solutions are stable at room temperature for two weeks. If a precipitate forms in Stop Solution during thawing, it should be redissolved by incubation at 37°C for 5 min. Do not shake the buffers vigorously.
 2. On the day of the assay first prepare the Standards.
 3. Bring vials with Buffer C and Assay Color Reagent to room temperature for thawing. It takes about 2-3 hours to melt. During this time perform extraction of the samples and prepare two aliquots of the sample extracts for separate measurement of NAD⁺ and NADH. Perform NAD⁺ and NADH assays on separate plates, one at the time.

EXTRACTION OF NAD⁺ AND NADH FROM BLOOD

1. Bring BUFFER A, NAD⁺ Stabilizing Reagent and NADH Stabilizing Reagent to room temperature before extraction.
2. Use fresh blood cooled on ice or thaw frozen blood samples in the ice-water bath for 12-15 min prior to the extraction. (**Tip:** use occasional quick cycles of warm hands during thawing time). We recommend extracting a maximum of eight samples at a time.
3. Pipette 500 μ L of BUFFER A into 1.5 mL test tubes.
4. Heat the aliquots of BUFFER A in a heating block for 1-2 min (max 5 min) at 75° - 80°C. Keep in the heating block until step 6.
5. Quickly mix blood samples by one up-and-down pipetting cycle before extraction, avoid foaming.
6. Add 100 μ L of blood by one quick firm move directly into the tube with hot BUFFER A kept in the heating block. Immediately mix by a few intensive up-and-down pipetting cycles with simultaneous rotation of the tip to efficiently mix the homogenate.
7. Incubate homogenate at 75° - 80°C for 1 min.
8. Cool the mixture in an ice-water bath for at least 5 min. After cooling on ice, homogenate should polymerize without any free liquid.
9. Positive Control does not contain proteins and does not require a hot extraction step. Thaw Positive Control for 10 min at room temperature. Mix 100 μ L of the Positive control with 500 μ L of non-heated BUFFER A. The obtained solution is ready for the next step of NAD⁺ and NADH stabilization.
10. Centrifuge at 20 000 x g for 10 min at 4°C. (Note: Positive control does not require centrifugation). Transfer the supernatant into a clean test tube and discard the pellet. Keep obtained total extracts containing both NAD⁺ and NADH at 4°C, covered with the foil lid until the next step. *Optionally*, the supernatants can be stored at -80°C for two weeks. Thaw on ice before preparation of aliquots in step 10.
11. Prepare two separate aliquots from the obtained extracts into clean test tubes - 150 μ L / tube.
12. To the first 150 μ L aliquot, add 100 μ L of **NAD⁺** Stabilizing Reagent to get Stabilized Extract with NAD⁺ (keeps NAD⁺ and removes NADH). Vortex and keep at room temperature. **Protect from light by covering with an aluminum foil lid.**
13. To the second 150 μ L aliquot, add 100 μ L of **NADH** Stabilizing Reagent to get Stabilized Extract with NADH (keeps NADH and removes NAD⁺). Vortex and incubate for 1-2 min at 75° - 80°C (max. 2 min). Cool on ice for 5 min, then keep at room temperature. **Protect from light by covering with an aluminum foil lid.**

NOTE: Final dilution of the initial whole blood sample will be 10 times.

NOTE: In case of supplementation with NAD precursors, the NAD⁺ Stabilized extract should be further diluted (two times) using MilliQ water before the assay. In this case, the dilution of the initial blood sample will be 20 times for NAD⁺. The NADH Stabilized extract does not require dilution.

PREPARATION OF STANDARDS

Handling notes: Protect the standards from light with aluminum foil.

1. Prepare standards on the day of the assay. Thaw tubes with 1mM standard stocks for 10 min on ice. Protect 1 mM stocks from light with a foil lid during thawing.
2. During thawing, prepare buffers for Standard dilutions:
 - **Buffer NAD⁺**: mix 3 mL of **BUFFER A** with 2 mL of **NAD⁺ stabilizing reagent**, vortex. Use to prepare **NAD⁺ Standards** in step 4.
 - **Buffer NADH**: mix 3 mL of **BUFFER A** with 2 mL of **NADH stabilizing reagent**, vortex. Use to prepare **NADH Standards** in step 5.
3. Prepare 50 μ M stock of NAD⁺ by adding 25 μ L of 1 mM NAD⁺ standard stock (provided) into 475 μ L of MilliQ water, vortex.
4. Prepare **NAD⁺ Standards** according to the scheme below by mixing indicated volumes of 50 μ M NAD⁺ stock and Buffer NAD⁺ for NAD⁺ Standards, vortex. Final volume of each Standard is 1mL.

STANDARD NUMBER	CONCENTRATION of NAD ⁺ (μ M)	STANDARD PREPARATION	
		50 M NAD ⁺ stock (μ L)	Buffer NAD ⁺ (μ L)
ST1/Plate Blank	0	0	1000
ST2	1	20	980
ST3	2	40	960
ST4	3	60	940
ST5	5	100	900

5. Prepare 10 μ M stock of NADH by adding 10 μ L of 1 mM NADH standard stock (provided) into 990 μ L of MilliQ water, vortex.
6. Prepare **NADH Standards** according to the scheme below by mixing indicated volumes of 10 μ M NADH stock and Buffer NADH for NADH Standards, vortex. Final volume of each Standard is 1 mL.

STANDARD NUMBER	CONCENTRATION of NADH (μ M)	STANDARD PREPARATION	
		10 M NADH stock (μ L)	Buffer NADH (μ L)
ST1/Plate Blank	0	0	1000
ST2	0.2	20	980
ST3	0.4	40	960
ST4	0.6	60	940
ST5	1	100	900

7. Cover the stand with ready Standards with an aluminum foil lid to protect from light and keep at 4°C prior to pipetting on the plate.

ASSAY PROCEDURE

Handling notes: Important! Protect the Master Mix from light (see practical considerations)

- Every assay contains Plate Blank and Sample Blanks to correct for all unspecific background signals. NAD⁺ or NADH Standards with concentration 0 μM (ST1) are used as Plate Blanks. Prepare separately four Sample Blank wells for the first four samples (BL UNK1-4 of Samples 1-4) to correct for unspecific interaction between the extract components and Assay Color reagent in Master Mix. Positive Control does not require a separate Blank.
 - ***All Sample Blanks are incubated with Master Mix WITHOUT added enzyme.***
 - Enzyme should be thawed just before addition into Master Mix. Brief centrifugation of the vial with enzyme at low speed prior to opening is recommended.
1. Pipette 20 μL of each NAD⁺ or NADH Standard in duplicates starting from ST1 (0 μM) according to the scheme on page 14.
 2. Pipette 20 μL of Positive Control and stabilized extracts in duplicates (see the scheme on p.14). For the first four samples, pipette one extra replicate to the indicated well (BL UNK1-4) according to the scheme on p.14. These four wells are Sample Blanks needed for the analysis without enzyme, to correct for unspecific interaction between the extract and Assay Color reagent within the Master Mix.
 3. Prepare Master Mix by adding Assay Color Reagent into BUFFER C, mix gently by rotation.
 4. Add 190 μL of Master Mix into each of the four Sample Blank wells.
 5. Add 40 μL of Enzyme into the remaining Master Mix. Mix gently, avoid foaming.
 6. Add 190 μL of Master Mix with added enzyme to all remaining wells including Positive Control, using a multichannel pipette, avoid foaming. Immediately cover the ready plate with the aluminum foil lid.
 7. For NAD⁺ assay: incubate for 4-6 min at room temperature.
For NADH assay: incubate for 8-10 min at room temperature.
 8. Stop the reactions by adding 10 μL of Stop Solution to each well in the same order as Master Mix using a multichannel pipette. Avoid foaming, gently shake the plate by hand on a table surface, remove bubbles.
 9. Measure light absorbance at 573 nm within 5-10 min after stopping the reactions. If possible, shake the plate inside the microplate reader for 5 sec before the measurement.

CALCULATION OF RESULTS

1. Read the plate at 573 nm within 5-10 min after stopping the reaction
2. Subtract the average of ST1 (0 μM NAD⁺ and 0 μM NADH) from the absorbance reading in each well on the plate.
3. Average the duplicate readings for each standard (ST1 – ST5). Create a standard curve by plotting the mean absorbance for each standard on the Y-axis against the concentration (in μM) on the X-axis and perform a simple linear regression fitting of the standard curve.
4. Average readings for four Sample Blanks (BL UNK1-4). The obtained value represents unspecific reactions of the extracts with Assay Color Reagent in the Master Mix.
5. Average the duplicate readings for each sample (UNK).
6. Average the duplicate readings for Positive control.
7. Find concentration of NAD⁺ and NADH from averaged absorbance readings of samples, Sample Blank and Positive control using formula of linear regression of the Standard curve. Correct obtained concentrations in Samples for Sample Blank and multiply on 10 to obtain the concentration (μM) of NAD⁺ or NADH metabolite in the initial blood sample. If NAD⁺ Stabilized extracts have been additionally diluted, the concentration must be multiplied by the dilution factor.

INTERNAL QUALITY CONTROL

Concentration of NAD⁺ in the positive control is expected to be within the range of 23 - 26 μM , concentration of NADH is expected to be within the range of 0.6 - 1.7 μM .


- If the results for the positive control are outside the range specified above, the results cannot be used unless a satisfactory explanation for the discrepancy has been given.
- Acceptance criteria for the difference between the duplicate results of the samples should rely on Good Laboratory Practices.
- It is a good practice to check visually the curve fit selected by the computer.

PERFORMANCE AND LIMITATIONS

A. TYPICAL DATA


The standard curve and the unknown concentration values are provided for demonstration only and should never be used instead of the real-time calibration curve.

STANDARD CURVE FOR NAD⁺

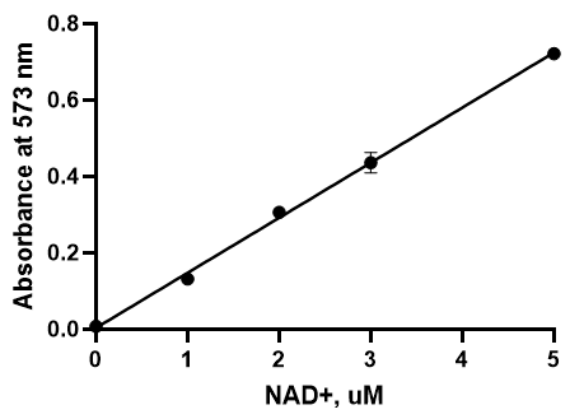

<i>Standard</i>	<i>NAD⁺ (μM)</i>	<i>Absorbance (573 nm)</i>	<i>Average Plate BL corrected</i>
ST1/Plate	0	0.051	0
Blank		0.051	
ST2	1	0.172 0.179	0.124
ST3	2	0.356 0.344	0.299
ST4	3	0.499 0.461	0.429
ST5	5	0.767 0.764	0.714

STANDARD CURVE FOR NADH

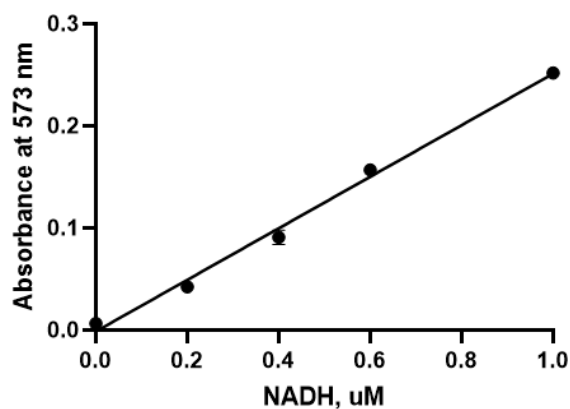

<i>Standard</i>	<i>NADH (μM)</i>	<i>Absorbance (573 nm)</i>	<i>Average Plate BL corrected</i>
ST1/Plate	0	0.060	0
Blank		0.059	
ST2	0.2	0.092 0.099	0.036
ST3	0.4	0.149 0.139	0.083
ST4	0.6	0.209 0.211	0.150
ST5	1	0.306 0.304	0.245

CALCULATION OF RESULTS FOR NAD⁺

Concentration values in the unknown samples (UNK) and Sample Blanks (BL UNK1-4) are determined from the linear fit formula of the NAD⁺ standard curve

Unknown	Concentration in Stabilized Extract (μM)	Average (μM)	Concentration in Stabilized Extract corrected by Sample Blank (BL UNK 1-4, μM)	Final NAD⁺ concentration in the original sample (μM)*
UNK 1	2.355	2.416	2.237	22.37
	2.428			
	2.464			
UNK 2	2.610	2.564	2.385	23.85
	2.592			
	2.492			
UNK 3	2.419	2.364	2.185	21.85
	2.355			
	2.319			
UNK 4	2.492	2.440	2.261	22.61
	2.391			
	2.437			
UNK 5	2.501	2.449	2.270	22.70
	2.446			
	2.401			
UNK 6	2.637	2.601	2.422	24.22
	2.610			
	2.555			
BL UNK 1	0.219	0.179	-	
BL UNK 2	0.163			
BL UNK 3	0.154			
BL UNK 4	0.179			

*Corrected by dilution factor x10

CALCULATION OF RESULTS FOR NADH

Concentration values in the unknown samples (UNK) and Sample Blanks (BL UNK1-4) are determined from the linear fit formula of the NADH standard curve.

Unknown	Concentration in Stabilized Extract (μM)	Average (μM)	Concentration in Stabilized Extract corrected by Sample Blank (BL UNK 1-4, μM)	Final NADH concentration in the original sample (μM)*
UNK 1	0.527	0.525	0.103	1.03
	0.527			
	0.522			
UNK 2	0.537	0.532	0.113	1.11
	0.532			
	0.527			
UNK 3	0.567	0.564	0.142	1.42
	0.562			
	0.562			
UNK 4	0.562	0.564	0.142	1.42
	0.562			
	0.567			
UNK 5	0.578	0.569	0.147	1.47
	0.562			
	0.567			
UNK 6	0.568	0.576	0.154	1.54
	0.578			
	0.583			
BL UNK 1	0.426	0.422	-	
BL UNK 2	0.410			
BL UNK 3	0.432			
BL UNK 4	0.422			

*Corrected by dilution factor x10

B. LIMITS OF DETECTION

The Limit of Blank (LOB) for Q-NADMED Blood is presented in the table below (LOB \pm Standard deviation (SD)).

Limit of Blank	
	pmol/well
NAD ⁺	1.84 \pm 0.9
NADH	2.10 \pm 0.5

The Limit of Detection (LOD) was calculated from NAD⁺ and NADH standard curves and is presented in the table below (LOD \pm SD).

Low detection limit	
	μ M in whole blood
NAD ⁺	0.33 \pm 0.2
NADH	0.19 \pm 0.05

The Limit of Quantitation (LOQ) is presented in the table below (LOQ \pm SD).

Limit of Quantitation	
	μ M in whole blood
NAD ⁺	0.66 \pm 0.3
NADH	0.40 \pm 0.1

A. PRECISION AND REPRODUCIBILITY

Precision of the assay performance was determined by Intra-assay variation in measurement. Intra-Assay precision is presented in the table below (CV = coefficient of variation).

Intra-assay precision	
	CV (%) \pm SD
NAD ⁺	1.48 \pm 0.8
NADH	3.33 \pm 1.5

Results of the assay reproducibility are summarized in the table below (N = number).

Reproducibility						
	NAD⁺			NADH		
Sample	Ctr1	Ctr2	Ctr3	Ctr1	Ctr2	Ctr3
N of measurements *	9	9	9	9	9	9
Mean (μ M)	27.41	29.41	22.00	0.55	0.71	0.64
Standard deviation	0.62	1.31	0.87	0.03	0.05	0.05
CV (%)	2.28	4.45	3.95	5.20	7.06	8.45

B. ACCURACY

The accuracy of the assay was calculated from samples with known amount of pure NAD⁺ and NADH. The results are summarized in the table below (assay accuracy \pm SD).

Accuracy (%)		
NAD ⁺	N = 32	97.13 \pm 7.6
NADH	N = 25	104.22 \pm 16.5

C. ASSAY CUT-OFF

The low and high cut-off values represent the smallest and the highest concentration observed in 5-7% individuals of given population extract. The cut-off values are summarized in the table below.


Cut-off value		
	Low	High
NAD⁺ (μ M)	20	36
NADH (μ M)	0.6	1.8

D. PERFORMANCE CHARACTERISTICS





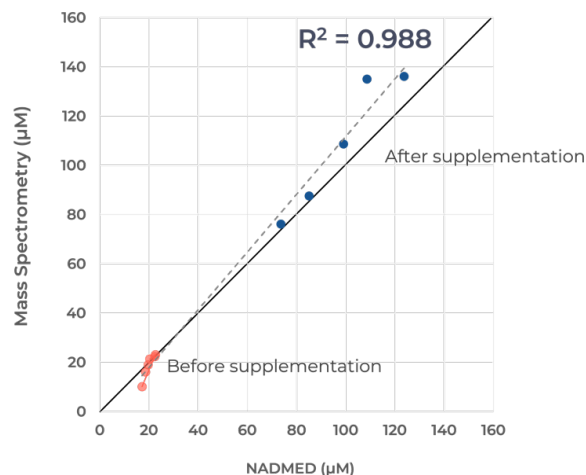


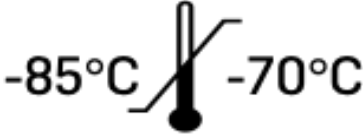
Interference of other metabolites in the extract was not separately investigated, as their contribution is low and taken into account by performing Blank analysis without added enzyme. **Warning:** presence of potassium sorbate, borate, pyridine, bismuth in a sample can cause enzyme inhibition, and thus causing underestimation of the results.






E. METHOD COMPARISON

To validate the performance of Q-NADMED, we measured NAD⁺ concentration in a set of control human blood samples that were also analyzed by Mass Spec. Frozen blood samples of five healthy subjects - before and after 16 weeks of niacin supplementation - were analyzed in parallel by Q-NADMED and Mass Spec. Results from Q-NADMED were concordant to those obtained by Mass Spec (Figure 10).




SYMBOLS

Symbol	
	Flammable liquid and vapor
	Warning/Danger
	Consult instructions for use
	Use-by date
	Catalogue number
	Batch code
	Manufacturer
	Upper limit of temperature

	<p>Do not use if package is damaged</p>
	<p>Keep dry</p>
	<p>Number or reactions</p>
	<p><i>In Vitro</i> Diagnostic Medical Device</p>
	<p>Protect from direct light</p>

SCHEMATIC PICTURES

EXTRACTION OF NAD+ AND NADH FROM BLOOD

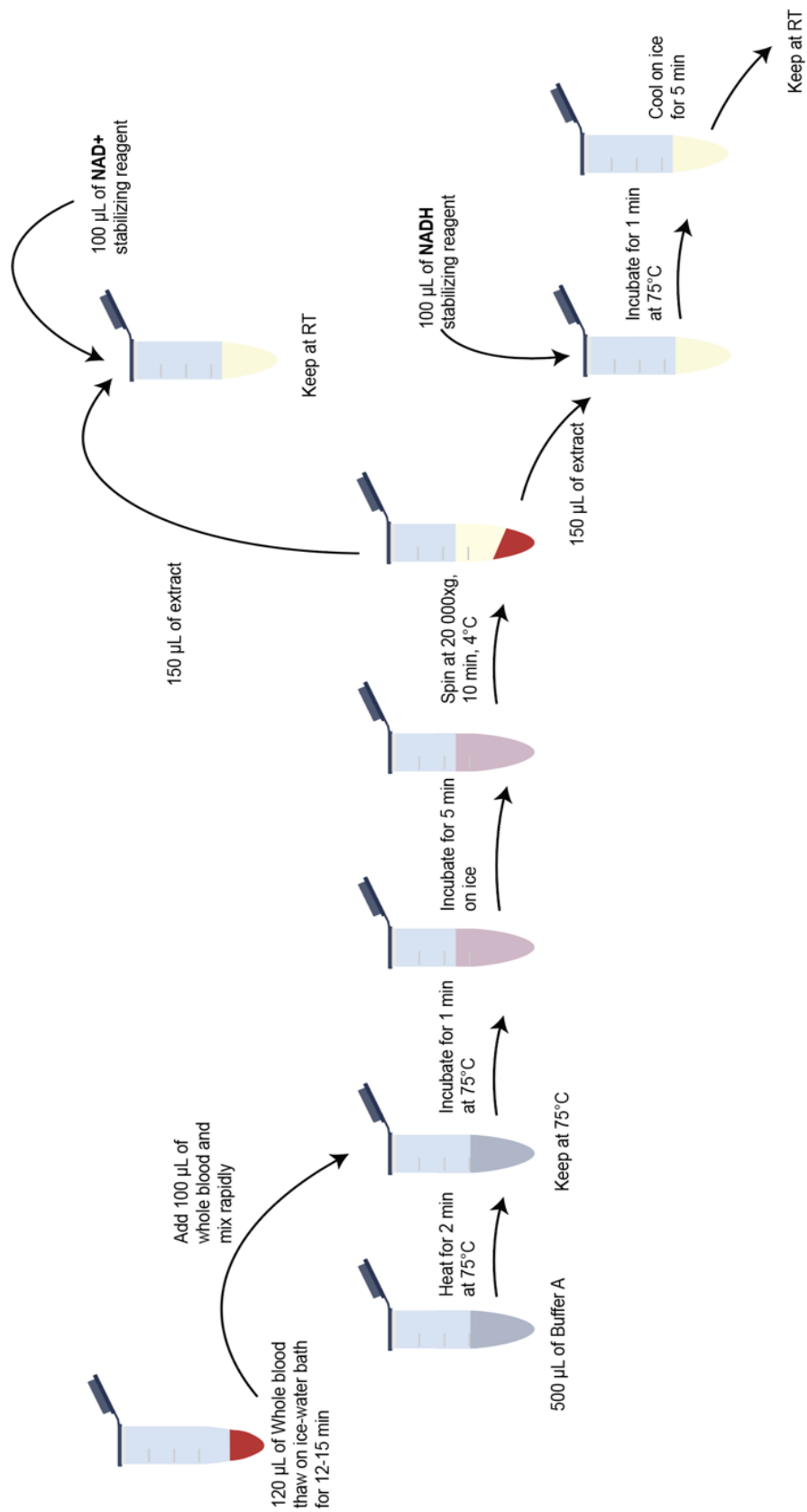

RECOMMENDED PLATE LAYOUT FOR NAD+ OR NADH MEASUREMENT

	1	2	3	4	5	6	7	8	9	10	11	12
A	St1 Plate BL	St1 Plate BL	UNK1	UNK1	UNK9	UNK9	UNK17	UNK17	UNK25	UNK25	UNK33	UNK33
B	St2	St2	UNK2	UNK2	UNK10	UNK10	UNK18	UNK18	UNK26	UNK26	UNK34	UNK34
C	St3	St3	UNK3	UNK3	UNK11	UNK11	UNK19	UNK19	UNK27	UNK27	UNK35	UNK35
D	St4	St4	UNK4	UNK4	UNK12	UNK12	UNK20	UNK20	UNK28	UNK28	UNK36	UNK36
E	St5	St5	UNK5	UNK5	UNK13	UNK13	UNK21	UNK21	UNK29	UNK29	UNK37	UNK37
F	PosCtr	PosCtr	UNK6	UNK6	UNK14	UNK14	UNK22	UNK22	UNK30	UNK30	UNK38	UNK38
G	BL UNK1	BL UNK2	UNK7	UNK7	UNK15	UNK15	UNK23	UNK23	UNK31	UNK31	UNK39	UNK39
H	BL UNK3	BL UNK4	UNK8	UNK8	UNK16	UNK16	UNK24	UNK24	UNK32	UNK32	UNK40	UNK40

Plate layout for NAD+ or NADH assay: St – standard, BL - blanks with indicated samples used for the blank reading, UNK – extracts from the samples with unknown metabolite concentration, PosCtr – Positive Control. Use 20 µL of Standards, Positive Control and Samples per well. Sample Blanks of the first four samples are analyzed in the Master Mix without added enzyme.

NOTES

PLATE LAYOUT

Use this plate layout to record standards and samples assayed.

12								
11								
10								
9								
8								
7								
6								
5								
4								
3								
2								
1								
	A	B	C	D	E	F	G	H